

STRAIGHT TO THE POINT WITH PATIENTS VISUAL OUTCOMES

Monofocal plus IOL




First IOL mimicking the optical
performance of the human lens

VOPTICA
SMART VISUAL OPTICS




ADVANCED OPTICAL TECHNOLOGY

The ArtIOLs® reproduce the imaging formation properties of the **natural crystalline lens** in the periphery of the retina with their **inverted meniscus design**¹

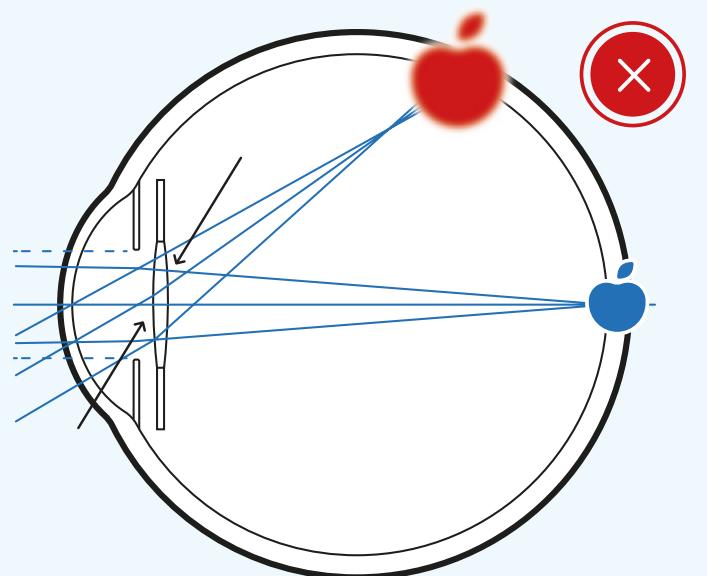


Standard IOL
(Biconvex)

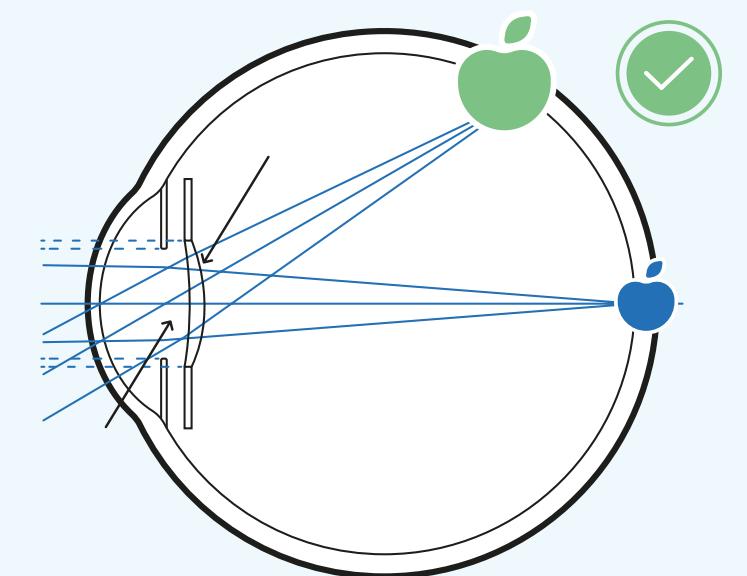


ArtIOLs®
(Inverted meniscus)

Patients implanted with ArtIOLs® present a **reduction of spherical equivalent by 1.97D at 40° of temporal eccentricity** compared to standard biconvex IOLs¹

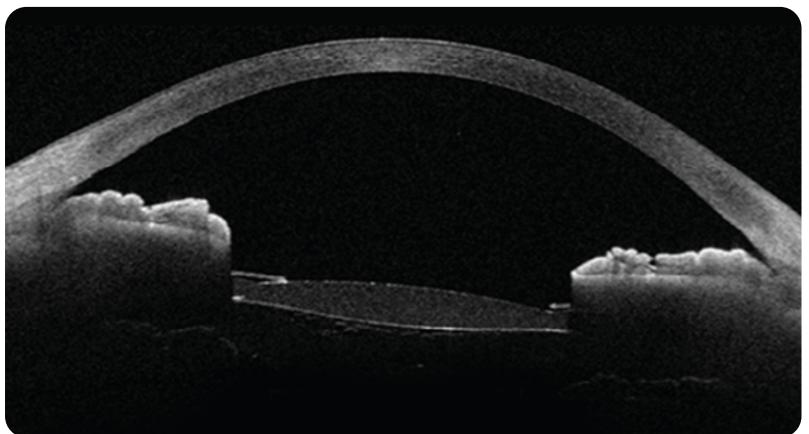

Biconvex IOL

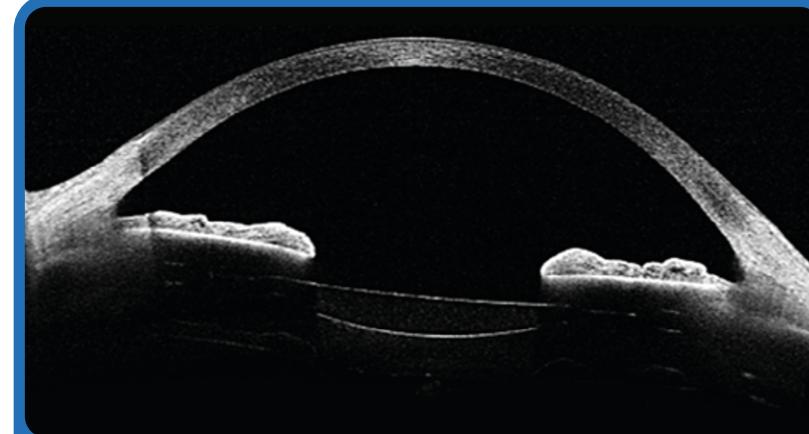

ArtIOLs®
Inverted Meniscus IOL

Anterior chamber OCT image of an implanted inverted meniscus ArtIOLs® in comparison with a standard biconvex IOL²



Standard IOL (Biconvex)⁴



ArtIOLs® (Inverted meniscus)

3. Harilaos G, et al. Quantification And Analysis Of Negative Dysphotopsia And Associated Phenomena In A Physical Model Of The Eye. Invest. Ophthalmol. Vis. Sci. 2022;63(7):1798 – F0414.

4. Pusnik A, et al. Dysphotopsias or Unwanted Visual Phenomena after Cataract Surgery. Life (Basel). 2022 Dec 24;13(1):53

1. Villegas EA et al. Peripheral, Refraction and Contrast Detection Sensitivity in Pseudophakic Patients Implanted With a New Meniscus Intraocular Lens. J Refract Surg. 2022 Apr;38(4):229-234

2. Marín JM, et al. Visual Performance at All Distances and Patient Satisfaction With a New Aspheric Inverted Meniscus Intraocular Lens. J Refract Surg. 2023 Sep;39(9):582-588.

PROPRIETARY ASPHERIC ANTERIOR AND POSTERIOR SURFACES

ArtIOLs provide a field of curvature to match the shape of the retina, achieving correction of myopic errors up to 2.00D at 40° of excentricity¹

Visual metrics

Central optics¹

Less peripheral defocus and astigmatism¹

Less peripheral distortion⁵

Superior Peripheral Contrast¹

Standard IOL (Biconvex)




ArtIOLs® (Inverted meniscus)



1. Villegas EA et al. Peripheral, Refraction and Contrast Detection Sensitivity in Pseudophakic Patients Implanted With a New Meniscus Intraocular Lens. J Refract Surg. 2022 Apr;38(4):229-234
5. Robles C, et al. Inverted meniscus IOLs reduce retinal distortion in the peripheral visual field. Invest. Ophthalmol. Vis. Sci. 2022;63(7):3508

ArtIOLs® 25 extend the depth of focus and provide superior visual outcomes at intermediate and near distances when compared to standard Monofocal IOL⁶

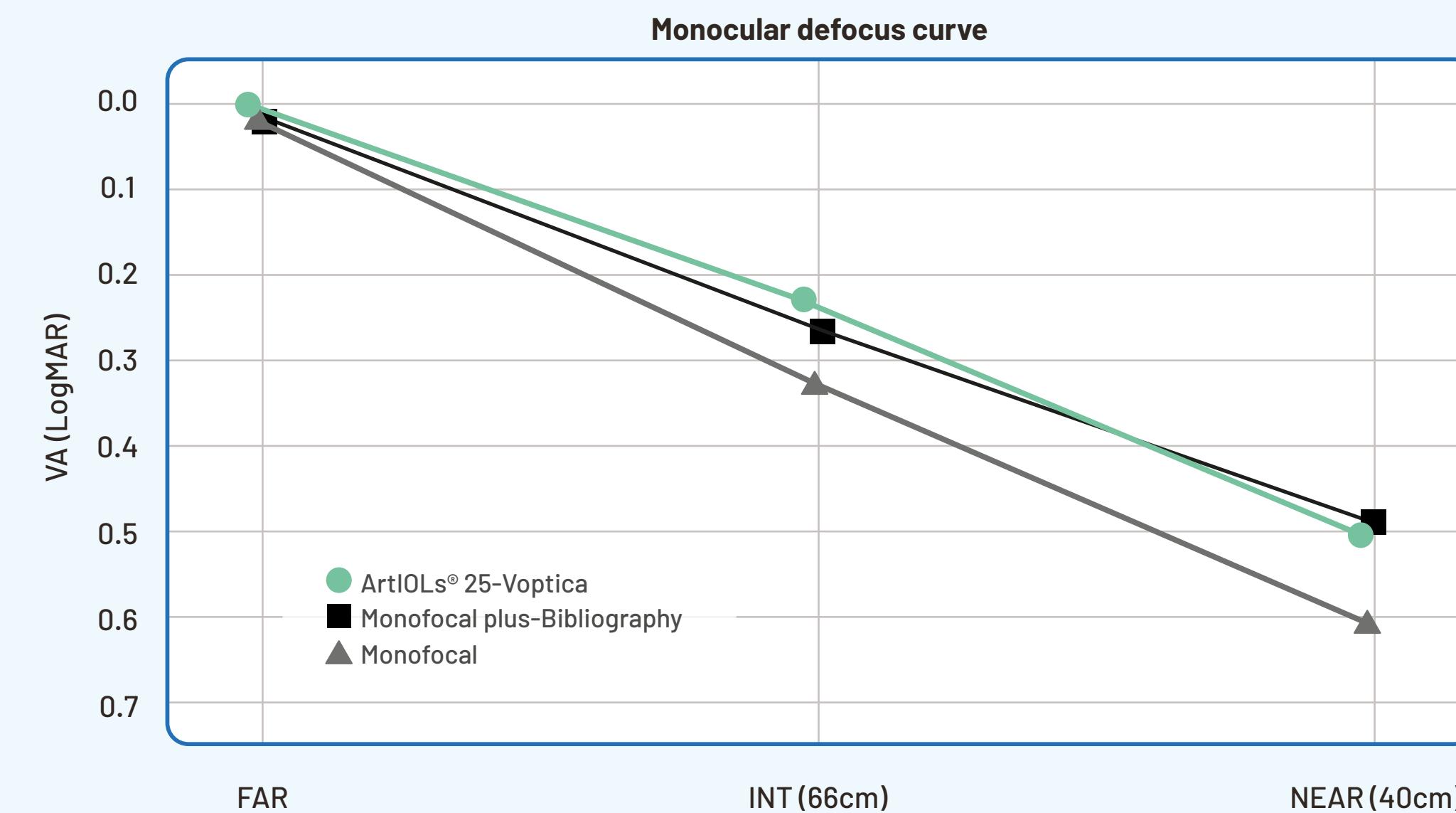

6. Voptica data on file: PRJ-INV_006

ArtIOLs® 25

TECHNICAL SPECIFICATIONS

voptica.com

Support:

info@voptica.com

(0034) 868 881 714



Optics

Lens type	Single piece foldable lens
Optical design	Aspheric optics with extended depth of focus
Shape	Inverted meniscus
Material	Hydrophobic Acrylic
	UV absorbing and blue light filter
Power range	+12.00 to +30.00 Diopters (0.50D steps)
Optical diameter	6.0 mm
Total diameter	13.0 mm
Refractive index	1.54 (glistening free)
Edge design	Square

Optical Biometry

Suggested A-Constant*	ArtIOLs® 25
SRK/T	120.0
Haigis	a0= 0.720, a1=0.297 and a2=0.169

Haptics

Haptic design	C-L haptic
---------------	------------

Delivery System

Injector-Cartridge set	Single use
Recommended incision size	≥ 2.2mm
Lens delivery	Single handed plunger

*It is recommended that surgeons personalize the constants they use.



2460
Coregroup Sight
Solution PL

MP-TR-013

VOPTICA
SMART VISUAL OPTICS